

Transportable compact wheelchair

Standard equipped with tension bar

Also available in paediatric sizes

Half folding backrest can be fitted as standard for UK customers.

Tension adjustable backrest.
Backrest height adjustable 42½, 45 and 47½ cm.

Seat cushion, 5 cm deep with non slip base (optional).

Removable footrests, swing away inwards/outwards. Easy to use grips for added comfort.

Tie down points illustrated for transportation.

Adjustable nylon heel loops on the footplates.

Height adjustable push handles.

Equipped with tension bar for added support.

Comfortable, flip backward and detachable armrests.

Multi-adjustable rear axle position.

Depth and angle adjustable braking system with extra comfy levers.

Quick release detachable rear wheel. S/P and transit. Lightweight spokes, fitted with high pressure performance tyres.

Castor adjustable to change seat height. Three options, 47½, 50, 52½ cm.



Multi-adjustable rear axle plate



Swing away and removable footrests.



Flip backwards and removable armrests.



Height adjustable footplates.

'Smart mobility solutions'


Standard specifications

Seat width: 40, 45, 50, 55 and 60 cm
 Seat depth: 43 and 46 cm
 Seat height: 47½, 50 and 52½ cm
 Backrest height: 42½, 45 and 47½ cm (adjustable)
 Total width: Seat width + 17 cm
 Total length: > 103 cm
 Total height: > 95 cm
 Width folded: 31 cm
 Armrest height: 26 or 19-35 cm adjustable (optional)

Front castor size: 7" x 1¼" (other sizes available)
 Rear wheel (sp): 24" x 1¾" (other sizes available)

Total weight (from): 16.9 kg (depending on the model)
 Weight (transport): 8.8 kg
 Max. user weight: 130 kg (seat width 50 and 55, 150 kg)

Frame colour available in:

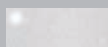

Product complies with crash test ISO 7176-19:2001 and EN12183



All rights reserved. Copyright 2008 Van Os Medical B.V.®

www.vanosmedical.com



Van Os Medical UK Ltd
 Corunna House
 42-44 Ousegate
 Selby
 YO8 4NH

Tel.: +44 (0) 1757 701177
 Fax: +44 (0) 1757 706011

sales@vanosmedical.com