

Lightweight foldable wheelchair

Available in eight seat widths
Standard equipped with quick release rear wheels

Tension adjustable seat and back upholstery.
Backrest height adjustable 42½, 45 and 47½ cm.

Armpads made of comfortable polyurethane,
contoured for extra comfort.
Flip backwards and removable armrests.

Seat cushion, 5 cm deep with
non slip base (optional).
Adjustable seat upholstery.

Footrests detachable,
swing away and swing in.
Height adjustable
footplates.

Tie down points
illustrated for
transportation.

Footplate angle
adjustable 90 - 105.
Adjustable nylon heel loops.

Height adjustable push handles.

Half folding backrest can be fitted as
standard for UK customers.
(For 55 and 60 cm seat width,
only fixed backrest available)

Depth and angle adjustable
braking system with extra
comfy levers.

Multi adjustable rear axle.
(Creates 3 options).

Quick release detachable rear
wheel. S/P and transit.
Lightweight spokes and fitted
with flat free PU tyre.
Equipped with aluminium
handrims.

Front castors 7" x 1 1/4" polyurethane
Multi position front fork, angle adjustable
front castor, horizontal adjustable.



Long armrests, height adjustable and removable (optional).



With:
- Seat height adjustable.
- Flat free (PU) front and rear wheels.
- Angle adjustable front fork.



Supremely comfortable length compensating
elevating legrest, includes easy adjustable calf
support (optional).



Flip-backward and detachable
armrests. Height adjustable with push
button (optional).

'Smart mobility solutions'


Standard specifications

Seat width: 35, 37½, 40, 42½, 45, 50, 55 and 60 cm
Seat depth: 43 and 46 cm
Seat height: 45, 47, 50, 52 and 54 cm
Backrest height: 42½, 45 and 47½ cm (adjustable)
Total width: Seat width + 17 cm
Total length: > 103 cm
Total height: > 95 cm
Width folded: 31 cm
Armrest height: 26 or 19-35 cm adjustable (optional)

Front castor size: 7" x 1 1/4" (other sizes available)
Rear wheel (sp): 24" x 1 3/8" (other sizes available)
Rear wheel (transit): 12 1/2" x 2 1/4"

Total weight (from): 13.9 kg (depending on the model)
Weight (transport): 8.4 kg
Max. user weight: 130 kg
(seat width 50 and 55, 150 kg
(seat width 60, 175 kg)

Frame colour
available in:

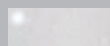

Product complies with crash test
ISO 7176-19:2001
and EN12183

All rights reserved. Copyright 2008® Van Os Medical B.V.®



Van Os Medical UK Ltd
Corunna House
42-44 Ousegate
Selby
YO8 4NH

Tel.: +44 (0) 1757 701177
Fax: +44 (0) 1757 706011

sales@vanosmedical.com