

Lightweight comfort wheelchair

Angle adjustable backrest and lightweight frame

Backrest height adjustable 42½, 45 and 47½ cm.
Foldable aluminum crossbrace.

Armpads made of comfortable polyurethane, contoured for extra comfort.
Height adjustable, flip backwards and removable armrests.

Seat cushion, 5 cm deep with non slip base.
Adjustable seat upholstery.

Quick release removable footrests.
Height adjustable footplates.

Tie down points illustrated for transportation.

Footplate angle adjustable 90 - 105.
Adjustable nylon heel loops.

Front castors 7" x 1 1/4" polyurethane
Multi position front fork, angle adjustable front castor, horizontal adjustable.

Height adjustable push handles, also angle adjustable.

Foldable tension bar, offering extra stability and greater comfort for the user.

Stepless adjustable backrest 0-45, operable with the control levers under the push handles.

Depth and angle adjustable aluminium brakes with extra comfy levers.

Quick release detachable rear wheels, S/P and transit.
Lightweight, fitted with PU tyre.

Anti-tip system, quickly height adjustable and detachable.
With integrated tipper system.



Also available in Transit version, fitted with 12½" rear wheels. Seat height standard is 52 cm.



With:
- Seat height adjustable.
- Flat free (PU) front and rear wheels.



Supremely comfortable length compensating elevating legrest, includes easy adjustable calf support (optional).



Equipped with stepless adjustable backrest 0-45.


Standard specifications

Seat width: 37½, 40, 42½, 45 and 50 cm
Seat depth: 43 cm
Seat height: 43, 45, 47, 50, 52 and 54 cm
Backrest height: 42½, 45 and 47½ cm (adjustable)
Total width: Seat width + 17 cm
Total length: > 103 cm
Total height: > 95 cm
Width folded: 31 cm
Armrest height: 22-28 cm adjustable (optional)

Front castor size: 7" x 1¼" (other sizes available)
Rear wheel (sp): 24" x 1¾" (other sizes available)
Rear wheel (transit): 12½" x 2¼"

Total weight (from): 14.9 kg (depending on the model)
Weight (transport): 9.4 kg
Max. user weight: 130 kg

Frame colour available in:

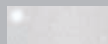

Product complies with crash test ISO 7176-19:2001 and EN12183



Van Os Medical UK Ltd
Corunna House
42-44 Ousegate
Selby
YO8 4NH

Tel.: +44 (0) 1757 701177
Fax: +44 (0) 1757 706011

sales@vanosmedical.com