

Excel

Navigator 3 'Compact'



All-round suspension compact mobility scooter

Equipped with an anti-tip system at the rear
Compact 3 wheel mobility scooter with extra comfort suspension

Comfortable head support for optimal use and comfort

Comfortable captains seat with:
An adjustable seat depth
Adjustable backrest
Lumbar adjustment

Folding armrest, width
angle and height adjustable

Removable collapsable
shopping basket
Including crutch holder

Detachable seat, four adjustable positions
for seat post
360 degrees rotation of seat

Easy to adjust rear suspension

Equipped with a 650 watt motor, equipped
with 4 motor brushes
Motor equipped with
electro-magnetic brake
Equipped with programmable digital
electronics, Dynamic Rhino 110 controller,
speed adjustable

Easy to use and understand dashboard
With large touch operation buttons
Setup so the user can use both their
left and right hand
Theft prevention system using ignition key
Speed adjustable in 6 steps

Indicators fitted, left and right flashing with
beep to inform of your intended movements

Height-adjustable steering column
Adjustable steering angle, stepless
adjustment, also can be folded
for transportation purposes

Impact-resistant plastic covers

Equipped with independent suspension
on the front and rear wheels



Equipped with an anti tip system.



The user can operate either left or right handed.



Equipped with front bumper for safety and to protect the plastic covers.



Easy to use waterproof dashboard, with touch operation.


'Smart mobility solutions'

Standard specifications

Overall width:	68 cm
Height rear:	120 cm
Height front:	114 cm
Overall length:	124 cm (including crutch holder)
Seat width:	47 cm
Seat height:	64-71 cm
Guardrail height:	20-30 cm
Total seat width:	58 cm
Wheelbase:	91 cm
Ground clearance:	16 cm
Speed:	Adjustable to 15 km/h (10 mph)*
Range:	Up to 50 km by 73 ah batteries

Total weight: 135 kg (depending on model)
Max. User Weight: 220 kg

Available in colours:

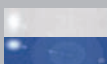

Product complies with crash test EN12184

We have a continuous product improvement policy and consequently reserve the right to amend design and specifications without prior notice. All errors and omissions excepted. All rights reserved. Copyright 2009 Van Os Medical B.V.



*For UK market scooters are only available with maximum speed 8 mph

www.vanosmedical.com



Van Os Medical UK Ltd
Corunna House
42-44 Ousegate
Selby
YO9 4NH

Tel.: +44 (0) 1757 701177
Fax: +44 (0) 1757 706011

sales@vanosmedical.com