

### Fixed frame sports wheelchair

Available in 3 seat widths

Aluminium lightweight sports wheelchair

Back upholstery on the back tubes.  
Padded back adjustable with tension straps.

Equipped with clothes protectors.

Aluminium 6 spoke wheels.  
Quick-release detachable rear wheels.  
Narrow tyres on the rear wheels.

Stepless front seat angle adjustment.  
Wheelbase adjustment through the ball head (2 positions possible).  
Quick-release detachable front forks.

Front castors 5" x 1" (solid).

Foldable back.  
Backrest angle adjustable.

Nylon calf support.



Super light, so extremely manoeuvrable

The ideal recreation sports wheelchair. The box frame guarantees an optimum stiffness and excellent driving characteristics.

Seat position is completely adjustable, the wheelbase, seat angle and seat height can be adapted to wishes of the user.



Standard equipped with angle adjustable aluminium brakes.



Standard equipped with design front forks solid 5" front wheels.


#### Standard specifications

Seat width: 40, 45 and 50 cm  
Seat depth: 43 cm  
Seat height: 48 and 50 cm  
Backrest height: 25 - 46 cm  
Total width: Seat width + 23 cm  
Total length: > 74 cm  
Total height: > 67 cm  
Seat angle: 6 - 17°  
Backrest angle: -5 till 15° stepless  
Wheel camber: 0, 3 and 6°

Front castor size: 5" x 1" (other sizes available)  
Rear wheel (sp): 24" x 1" (other sizes available)

Total weight (from): 10.9 kg (depending on the model)  
Weight (transport): 7 kg  
Max. user weight: 130 kg (seat width 50 cm, 150 kg)

Frame colour available in:

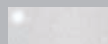

Product complies with crash test EN12183

All rights reserved. Copyright 2008© Van Os Medical B.V.®

[www.vanosmedical.com](http://www.vanosmedical.com)



Van Os Medical UK Ltd  
Corunna House  
42-44 Ousegate  
Selby  
YO8 4NH

Tel.: +44 (0) 1757 701177  
Fax: +44 (0) 1757 706011

[sales@vanosmedical.com](mailto:sales@vanosmedical.com)